

Euronda Pro System

Product, Process, Protection

EN



Euronda[®]
Pro System

Euronda®

It's all about protection

Founded in 1975, Euronda is now an Italian company with a global outlook, synonymous with outstanding protection in the medical field. With more than 200 employees, a 15,000 sq.m. headquarters, 4 subsidiaries abroad and affiliates worldwide, it is present on the market with 4 product lines dedicated to the safety of the dental environment and the protection of professionals, dentists and patients.

Monoart®, the range of disposable products that teams maximum quality with vibrant colours. **Pro System**, the comprehensive sterilization system that optimises and tracks the entire process. **Alle®**, the disposable line to create a sterile surgical field and prevent infection. **Eurosteril®**, rolls and pouches that ensure sterility is maintained over time.

Driven by a dynamic team, Euronda combines in-depth knowledge of the market with ongoing innovation in materials, processes and technologies to meet the needs of the sector and optimize work phases in full compliance with regulations, with particular attention to the environmental impact of the products and processes.





Accessorized

Equipped with stainless steel perforated basket and anti-deformation cover with non-drip system as standard. Optional beaker and beaker holder to wash small instruments such as burs or prostheses.



Clear

Intuitive and user-friendly thanks to the LED-backlit controller and light path that provides information on the status of the wash, guiding operators at a glance, even from a distance.



Ergonomic

Practical side handles and rear tap for easy movement and emptying.



Versatile

Users can choose from 6 pre-set cycles, or customise the cycle by setting the temperature up to 60°C and the time from 10 to 40 minutes.

Euronda
Pro System



Eurosonic® 4D

The new dimension of ultrasonic washing



Eurosonic® 4D is the most advanced ultrasonic tank in the Euronada Pro System range. Equipped with an LED-backlit display controller, it is intuitive and user-friendly thanks to the lights indicating the different cycles, times and temperatures. Durable and reliable thanks to its stainless steel tank (AISI 304L), Eurosonic® 4D allows a choice of 6 pre-set cycles and a free cycle that can be customised when it comes to time and temperature. The pre-set cycles work through 5 phases: heating of the liquid, degassing, removal of coarser particles, cleaning of the material through the combined action of micro-brushing and detergents, and a final phase that evenly distributes the ultrasonic waves for more uniform cleaning.



Size

33.5 x 26 x 20.5 cm (wxhxd)



Weight

5 Kg (with basket and cover)



Capacity

3.8 L tank



Supply voltage

220-230 V | 50/60 Hz

Power consumption

300 W



Pre-set cycles

- Instruments: 20 min at 40°C
- Burs: 10 min at 40°C
- Impression trays: 40 min at 60°C
- Cement removal: 10 min at 30°C
- Plaster removal: 30 min at 60°C
- Prostheses: 40 min at 60°C

Accessories included



Anti-deformation non-drip tinted cover



Perforated stainless steel basket

Optional accessories



Beaker holder 3 x 0.6 L



Beaker holder 1 x 0.6 L



0.6 L glass beaker with plastic lid



Simple

Eurosonic® 3D allows manual adjustment of the cleaning cycle time from 5 to 30 minutes. Users can also set the working temperature choosing between 60°C (default) and 35°C or start a cycle without heating.



Accessorized

The sound-absorbing cover with non-drip system protects operators from possible splashes and aerosols, preventing contamination of the work area and attenuating the noise produced by ultrasonic waves.



Smart

The central LED provides immediate information on the status of the cleaning cycle at any time.



Effective

The Sweep Mode technology increases the mechanical effectiveness of cavitation and ensures a more uniform cleaning due to a homogeneous distribution of ultrasonic waves within the tank.



Eurosonic® 3D

The most balanced choice for ultrasonic washing



Ultrasonic washing tank with pre-set working temperature of 60°C and adjustable time up to 30 minutes. Eurosonic® 3D offers maximum user-friendliness thanks to practical front buttons and luminous feedback on cycle status and correct operation. Its sturdy stainless steel tank (AISI 304L) and convenient side handles facilitate emptying, while the Sweep Mode technology ensures an effective and even cleaning.



Size

27 x 24 x 19.5 cm (wxhxd)



Weight

3.7 Kg (with basket and cover)



Capacity

2.7 L tank




Supply voltage

220-230 V | 50/60 Hz

Power consumption

280 W

Accessories included

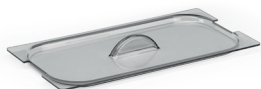

Anti-deformation non-drip tinted cover



Perforated stainless steel basket



0.6 L glass beaker with plastic lid



Beaker holder 2 x 0.6 L

Eurosonic[®] Energy

Safe and accessible ultrasonic washing



Ultrasonic tank with mechanical timer that allows users to set the duration of the cleaning cycle from 5 to 30 minutes. Eurosonic[®] Energy is equipped with the Sweep Mode technology, which evenly distributes the ultrasonic waves for deep and uniform cleaning. Sturdy with a stainless steel body (AISI 304L) and an anti-deformation non-drip cover, it has side handles for easy emptying and an effective and intuitive control panel.

Accessories included

- Perforated stainless steel basket;
- Sound-absorbing non-drip cover.

Optional accessories

- 0.6 L beaker with plastic lid;
- Beaker holder 2 x 0.6 L.



Size

27 x 24 x 22 cm (wxhxd)



Weight

3.7 Kg (with basket and cover)



Capacity

2.7 L tank



Supply voltage

230 V / 50/60 Hz

Power consumption

80 W



Speed

With Eurosonic[®] Energy, users can start cleaning instruments without the need to preheat the tank.



Smart

Indicator lights signal the status of the cleaning cycle or any errors and anomalies in its use at any time.



Simple

A single gesture is enough to start the cleaning cycle quickly and intuitively.

Euronda Spa® 2023

This catalogue may not be copied, reproduced or scanned by any means. All rights reserved.

Images, descriptions, quotations, dimensions and technical characteristics are to be considered purely indicative.

While maintaining the basic characteristics of the product, the manufacturer reserves the right to make changes at any time for the sole purpose of improvement.

Euronda Spa

Via Chizzalunga 1
36066 Sandrigo
Vicenza, Italy
t +39 0444 656111
f +39 0444 656199
m info@euronda.com
www.euronda.it

Deutschland

Euronda Deutschland
Am Landwehrbach 5
48341 Altenberge, Deutschland
t +49 2505 9389 0
f +49 2505 9389 29
m info@euronda.de
www.euronda.de

France

Euronda France Eurl
ZAC les Vallées, Avenue de Bruxelles
60110 Amblainville, France
t + 33 (0) 3 44 06 69 70
f + 33 (0) 3 44 02 03 89
m info@euronda.fr
www.euronda.fr

España

Identyd, Calle Músico
Antonio Rodríguez de Hita, 12
Bajo 30007 Murcia
t 900102034
m info@euronda.es
www.euronda.es

Russia

OOO Euronda Russia
193149 Leningradsakaya obl.,
Vsevolozhsky r-on,
d. Novosaratovka, ul.
Pokrovskaya, dom 41, litera B. Russia
t + 7 (812) 635 88 94
m info@eurondarussia.com
www.eurondarussia.com

Euronda[®]