

# Euronda Pro System

Product, Process,  
Protection

EN



**Euronda**<sup>®</sup>  
Pro System

# Euronda®

## It's all about protection

Founded in 1975, Euronda is now an Italian company with a global outlook, synonymous with outstanding protection in the medical field. With more than 200 employees, a 15,000 sq.m. headquarters, 4 subsidiaries abroad and affiliates worldwide, it is present on the market with 4 product lines dedicated to the safety of the dental environment and the protection of professionals, dentists and patients.

**Monoart®**, the range of disposable products that teams maximum quality with vibrant colours. **Pro System**, the comprehensive sterilization system that optimises and tracks the entire process. **Alle®**, the disposable line to create a sterile surgical field and prevent infection. **Eurosteril®**, rolls and pouches that ensure sterility is maintained over time.

Driven by a dynamic team, Euronda combines in-depth knowledge of the market with ongoing innovation in materials, processes and technologies to meet the needs of the sector and optimize work phases in full compliance with regulations, with particular attention to the environmental impact of the products and processes.





## Adjustable

The temperature can be adjusted from the control panel to ensure the perfect sealing with any type of material.



## Safe

An audible alarm alerts operators to any matters requiring their attention and the stand-by function automatically activates when the machine is not in use.




## Traceable

The built-in printer prints data directly onto the pouches. All information is transferable to a PC via USB, SD card and standard Ethernet connection thanks to the free E-Memory software.



## Accessories available

A comprehensive range of accessories are available to enhance usability, including a front sliding table and a table-top roll-holder with cutter.

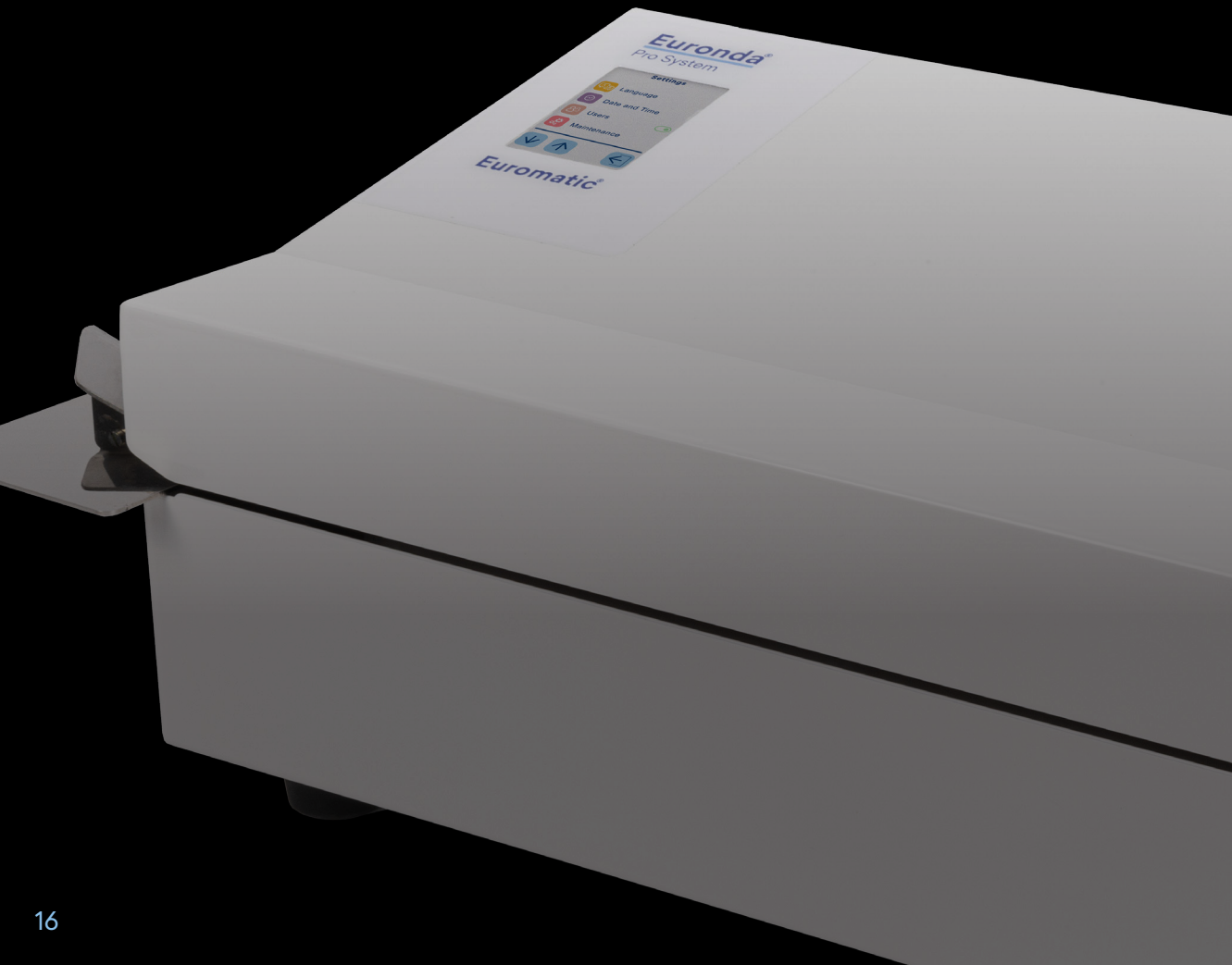

# Euromatic®

Simplifies, speeds up and keeps track of work



Euromatic® is the ideal thermosealing machine for practices with large quantities of pouches to seal. It validates every sealed pouch, ensuring maximum time savings due to multi-band sealing at a speed of 8 metres per minute. The large 2.8" touch screen display features a user-friendly interface that guides users through the menu options quickly and clearly. In addition, the screen undergoes a special treatment that repels and inhibits the proliferation of bacteria. The built-in printer prints useful parameters for traceability (production data, batch, operator, expiry date) directly on the pouch. All the sealing information can be transferred to a PC via a USB device, SD card and Ethernet connection as standard.



## Size

51 x 21.5 x 24 cm (wxhxd)



## Weight

14 Kg



## Sealing width

1.25 cm (double band)



## Supply voltage

200-240 V | 50/60 Hz

## Power consumption

500 W



## Sealing temperature


130 - 220°C




## Connectivity

Included: Ethernet, SD card, USB and E-Memory software


## Optional accessories

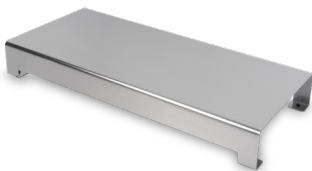

Height-adjustable front slide  
48 x 6 x 24 cm (wxhxd)

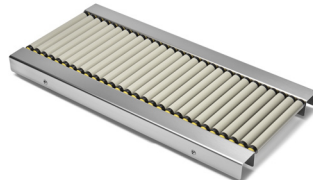

Height-adjustable front roller table  
79 x 7 x 29 cm (wxhxd)

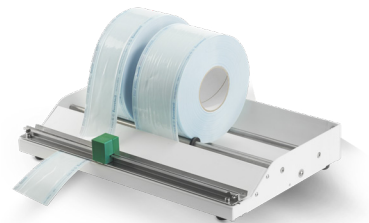

Bench-top roll holder and cutter unit  
48 x 8 x 30 cm (wxhxd)



## User-Friendly

The display makes it easy to choose the sealing temperature and provides real-time information on the effectiveness of operations. Thanks to sound signals, it minimises the risk of errors and accidents.



## Validating

It validates the sealing, confirming its success and storing the relevant data (time, temperature, pressure) on an SD card. The data can then be transferred to a PC via Ethernet or WiFi (optional).




## Practical

Ergonomic support surface and polycarbonate cover for constant visual inspection of the sealing area.



## Versatile

The patented double-level, height-adjustable roll holder can be fixed to a wall, minimising the space required for its use

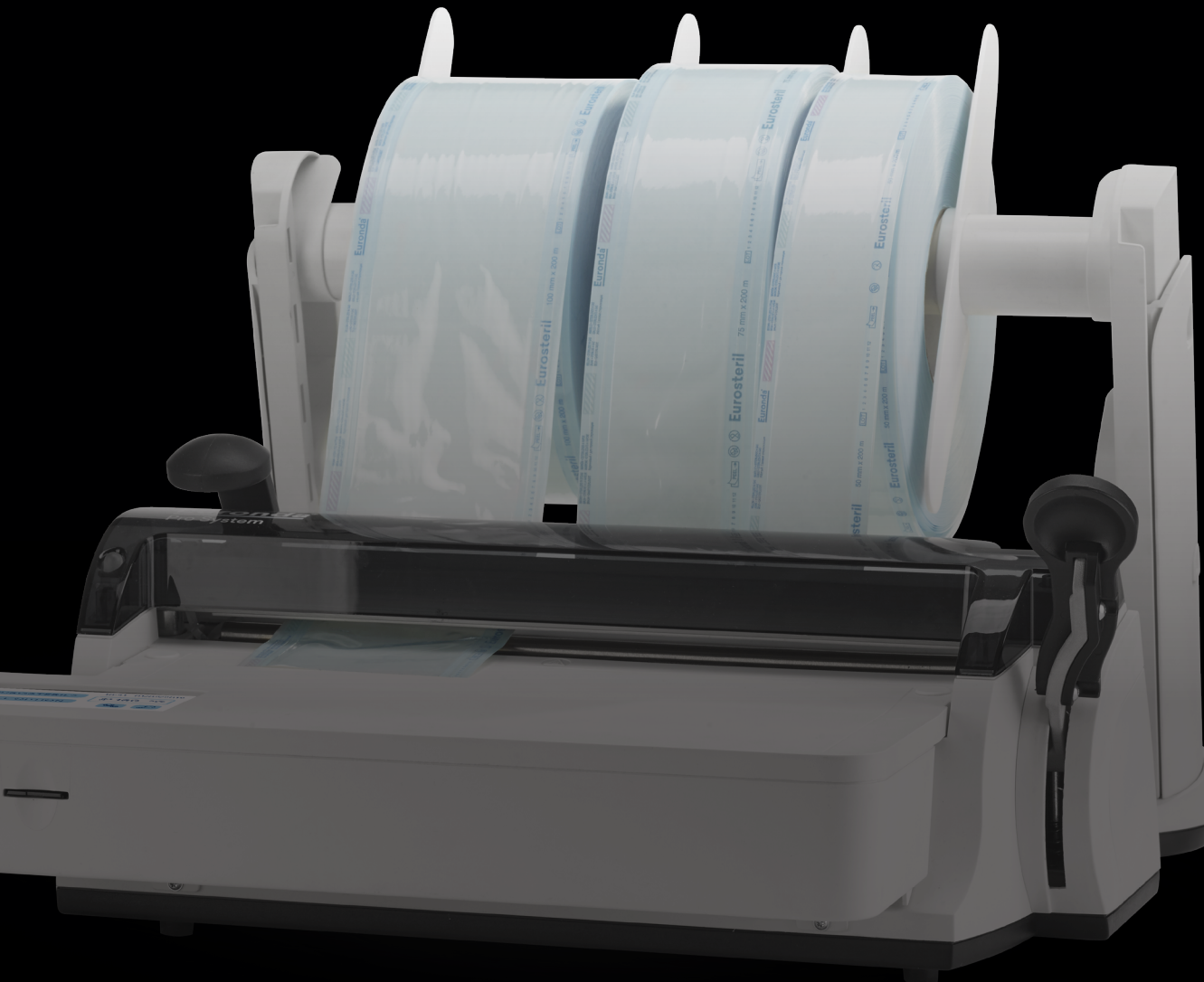

# Euroseal® Valida

## Digital traceability at every sealing



Euroseal® Valida can validate each sealing, saving the data concerning the cycle and the results on an SD card. It can be connected directly to a PC using Ethernet or WiFi (optional) connections for an immediate transfer of data. With its ergonomic design, it guarantees sturdiness and versatility, but also safety and ease of use. The built-in cutter and sealing lever located on the front optimize ergonomics and safety and make packing and sealing operations easier. The user experience is enhanced by the 4" colour touch screen display with user-friendly interface that allows users to quickly and intuitively set the language, date, time, operators, as well as the sealing temperature.



### Size (wxhxd)

- Standard: 50 x 31.5 x 49 cm
- With wall-mounted roll holder: 50 x 20 x 49 mm
- With double roll holder: 50 x 42 x 49 cm



### Weight

8.3 Kg



### Sealing width

1.2 cm



### Supply voltage

200-240 V | 50/60 Hz

### Power consumption

180 W



### Sealing temperature

170 - 200°C



### Connectivity


Included: Ethernet, SD card and E-Memory software  
Optional: WiFi



4" colour touch screen display



Improved position of Ethernet port

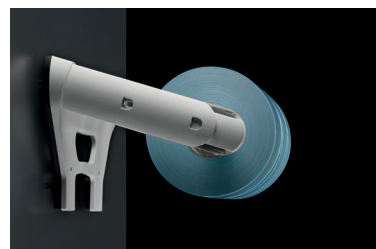

Roll holders can be wall mounted to optimise space



## Versatile

Patented height-adjustable roll holder that can be wall-mounted.



## Sturdy

Lightweight cutting-edge materials, built to withstand heavy usage.




## Practical

Spacious support surface that makes work practical, comfortable and easy.



## Safe

Anti roll-back blocking and scald prevention system, combined with sound and light signals.

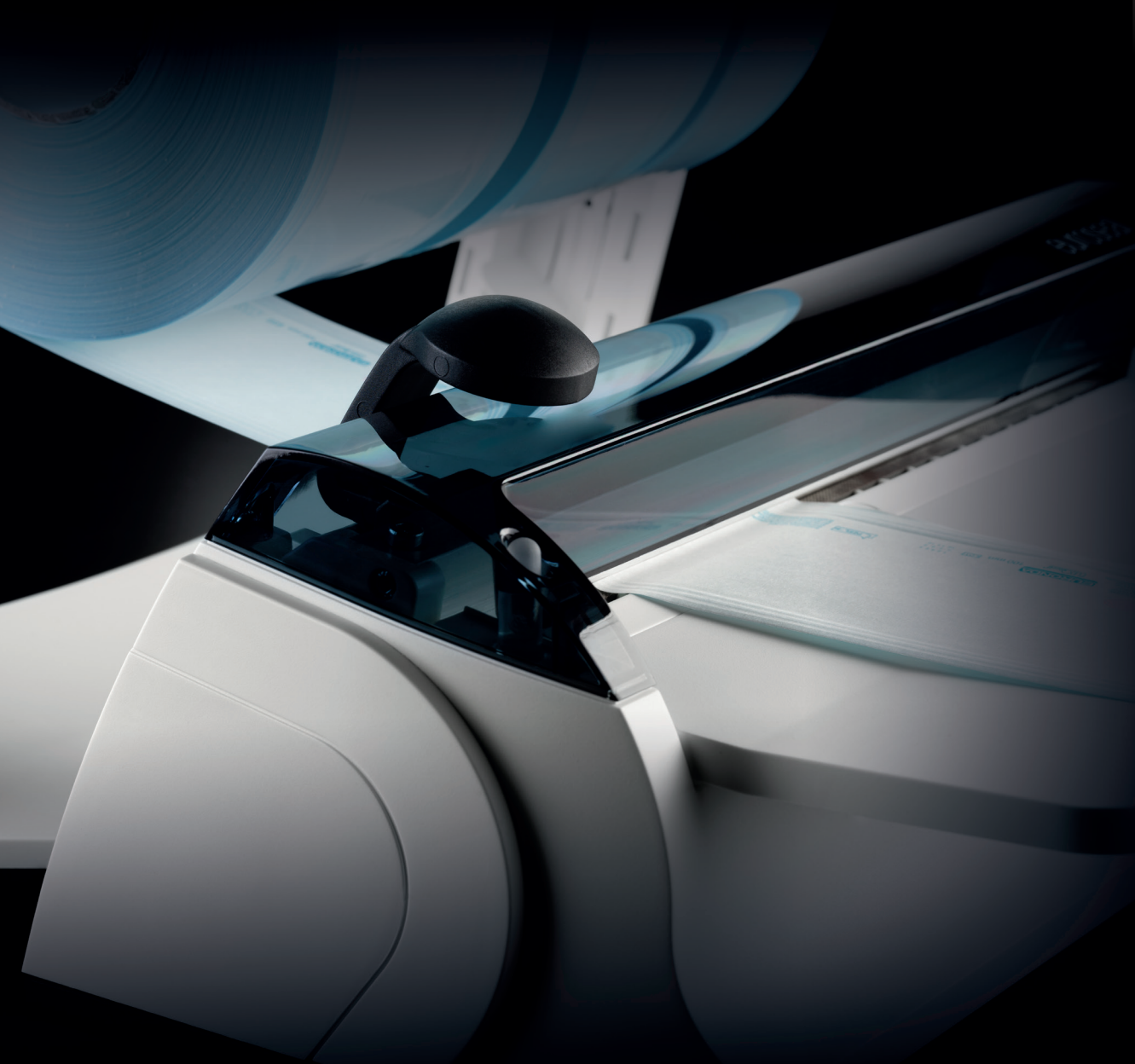

# Euroseal®

## The user-friendly sealer



Euroseal® is attractive, ergonomic and easy to use. The see-through polycarbonate cover allows operators to perform sealing operations with greater precision. Its special patented roll holder is adjustable in height and depth, and can also be mounted on the wall to suit all work spaces. The integrated support surface aids the sealing process and provides more practical working conditions. Anti roll-back blocking system.



### Size (wxhxd)

- Standard: 46 x 31.5 x 39 cm
- With wall-mounted roll holder: 46 x 20 x 30 mm
- With double roll holder: 46 x 42 x 39 cm



### Weight

7 Kg



### Sealing width

1.2 cm



### Supply voltage

200-240 V | 50/60 Hz

### Power consumption

100 W



### Sealing temperature

170 - 200°C



Integrated support surface



Adjustable height and depth of roll holder



## Ergonomic

The front cutter and sealing lever optimize ergonomics and safety, making packing and sealing operations easier.



## Compact

The machine is compact so it can be comfortably placed on a table top or vertically thanks to a kit for wall installation.



## Versatile

It enables 2 roll holder positions: high for best performance when the machine is wall-mounted, or low to speed up roll selection and changing operations.




## Sturdy

The quality of the materials and sturdy design make it practically indestructible, guaranteeing long-lasting efficiency and reliability.



# Euroseal® Infinity

A sealing machine  
is for forever

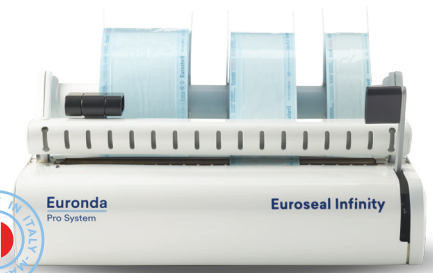

Euroseal® Infinity is practical, sturdy, reliable and easy to use. Its compact and functional design has been studied in detail to minimise the space required for its use; a feature implemented by the possibility of mounting the sealer to the wall for vertical use. It features a system that always maintains the roll inside the sealing area. The smooth surface facilitates and improves cleaning.



### Size

43 x 21.5 x 43 cm (wxhxd)



### Weight

7 Kg



### Sealing width

1.2 cm



### Supply voltage

200-240 V | 50/60 Hz

### Power consumption

180 W



### Sealing temperature

170 - 200°C



Easy-to-change integrated cutter



The roll holder can be positioned at 2 heights



Optional kit for wall installation

## **Euronda Spa® 2023**

This catalogue may not be copied, reproduced or scanned by any means. All rights reserved.

Images, descriptions, quotations, dimensions and technical characteristics are to be considered purely indicative.

While maintaining the basic characteristics of the product, the manufacturer reserves the right to make changes at any time for the sole purpose of improvement.

**Euronda Spa**

Via Chizzalunga 1  
36066 Sandrigo  
Vicenza, Italy  
**t** +39 0444 656111  
**f** +39 0444 656199  
**m** info@euronda.com  
**www.euronda.it**

**Deutschland**

Euronda Deutschland  
Am Landwehrbach 5  
48341 Altenberge, Deutschland  
**t** +49 2505 9389 0  
**f** +49 2505 9389 29  
**m** info@euronda.de  
**www.euronda.de**

**France**

Euronda France Eurl  
ZAC les Vallées, Avenue de Bruxelles  
60110 Amblainville, France  
**t** + 33 (0) 3 44 06 69 70  
**f** + 33 (0) 3 44 02 03 89  
**m** info@euronda.fr  
**www.euronda.fr**

**España**

Identyd, Calle Músico  
Antonio Rodríguez de Hita, 12  
Bajo 30007 Murcia  
**t** 900102034  
**m** info@euronda.es  
**www.euronda.es**

**Russia**

OOO Euronda Russia  
193149 Leningradsakaya obl.,  
Vsevolozhsky r-on,  
d. Novosaratovka, ul.  
Pokrovskaya, dom 41, litera B. Russia  
**t** + 7 (812) 635 88 94  
**m** info@eurondarussia.com  
**www.eurondarussia.com**

# **Euronda**<sup>®</sup>